

TIARIS™ Unstained Protein Marker

(Protein Research)

Ordering info

TBR0295, 2 x 250 µL

Description

TIARIS™ Unstained Protein Marker is a ready to use molecular weight marker. It is composed by 9 proteins ranging from 12 to 120 kDa. 12 and 50 kDa bands have increased intensity relative to the other bands for easy identification. The product includes a 12 kDa prestained band for electrophoresis and transfer monitoring.

Features

- Ready to use.
- Intense and accurate.
- Included highlighted band as visual reference.

Applications

- Size determination of proteins on SDS-PAGE.


Storage

Store at -20°C.

The product is shipped at room temperature (stable 3 months).

Quality Control

Band location and intensity by SDS-PAGE gel electrophoresis.

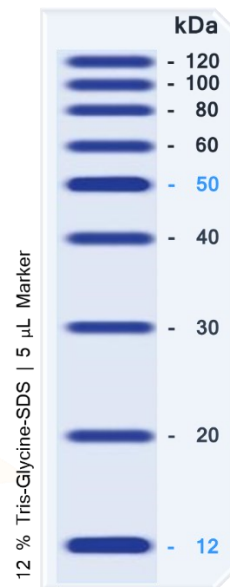

Use 3-5 µL of Tiaris™ Unstained Protein Marker per well.

Also available:

- Superfast-Blue™ Coomassie Stain (TBB0419)
- Acrylamide/ Bisacrylamide Solution (30%, 37.5 : 1) (TBR0135, TBR0136)
- Tris-Glycine-SDS Buffer 1x (TBB0335, TBB0336)
- Tris-Tricine SDS Buffer 10x (TBB0395, TBB0396)
- Ammonium Persulfate (TBR0140, TBR041)

TIARIS™ Unstained Protein Marker

(Investigación en Proteínas)

Referencias

TBR0295, 2 x 250 µL

Descripción

TIARIS™ Unstained Protein Marker is un marcador de peso molecular listo para usar. Está compuesto por 9 proteínas que van desde 12 hasta 120 kDa. Las bandas de 12 y 50 kDa tienen una mayor intensidad en comparación con las otras bandas para facilitar su identificación. El producto incluye una banda de 12 kDa preteñida para el monitoreo de la electroforesis y la transferencia.

Características

- Listo para usar.
- Intenso y preciso.
- Incluye bandas más intensas como referencia visual.

Aplicaciones

- Determinación del tamaño de proteínas en geles de SDS-PAGE.


Almacenaje

Conservar a -20°C.

El producto es enviado a temperatura ambiente (estable durante 3 meses).

Control de Calidad

Migración electroforética e intensidad de las bandas en gel SDS-PAGE.

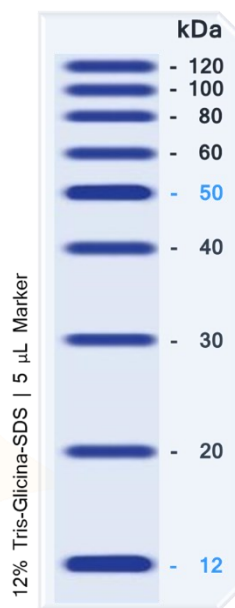

Usar 3-5 µL of Tiaris™ Unstained Protein Marker por pocillo.

También disponemos de:

- Superfast-Blue™ Coomassie Stain (TBB0419)
- Acrylamide/ Bisacrylamide Solution (30%, 37,5 : 1) (TBR0135, TBR0136)
- Tris-Glycine-SDS Buffer 1x (TBB0335, TBB0336)
- Tris-Tricine SDS Buffer 10x (TBB0395, TBB0396)
- Ammonium Persulfate (TBR0140, TBR041)